



Programming Crest Elite Digital™ Pillow Speakers

1. Locate your brand of TV in the **Entertainment Device Codes** section located on page two. You will be using the three-digit codes listed for your brand in step six of the programming procedure, if your pillow speaker has 10-key overlay. Four-digit device codes will be used for big button overlay. If your TV manufacturer isn't listed, please use code for the **AutoSense™**.

Note:

Please review the **Big Button Overlay Template** located to the right to decipher which keys correlate to the numbers 0, 1, 2 and 3. A four digit code will be required for this overlay as shown in the **Entertainment Device Code** section.

2. Plug the pillow speaker into the station.

3. Press and hold the **CC** and **mute** buttons.

4. Release the **mute** button after five seconds. Continue to hold the **CC** button.

5. Release the **CC** button after five seconds. You will have roughly 15 seconds to begin entering the three-digit device code (four-digit code if using the Big Button Overlay).


6. Enter the device code.

7. Test all the buttons on the pillow speaker. If the code doesn't work, program the pillow speaker by repeating steps two through five with another device code listed for your TV brand.

8. If none of the codes in the list work for your TV, please contact our Technical Support Department at **1-800-328-8908**.

Note:

The remote will return to normal operation throughout steps two through four if you simply release both the **CC** and **mute** keys. The remote will default back to the codes that were being used prior to programming.

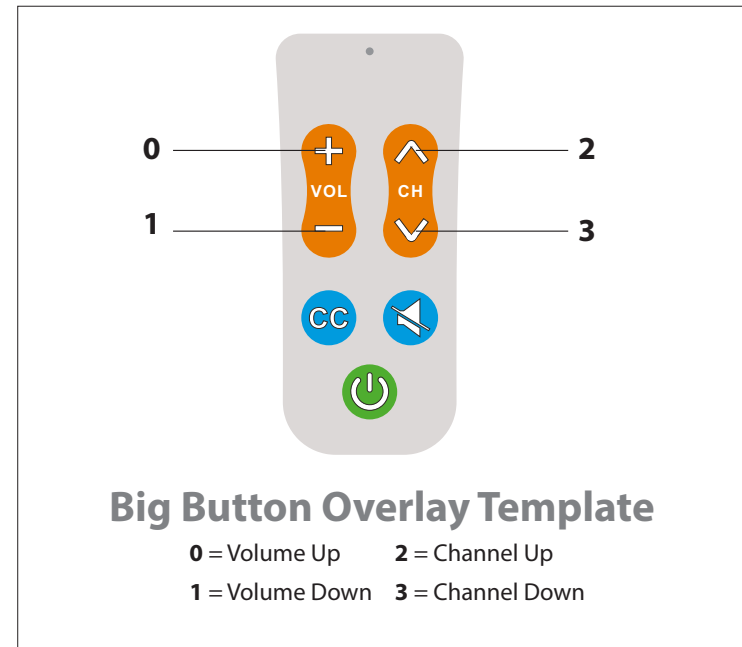

Entertainment Device Codes

Listed below are the most common TV brands used in healthcare facilities. If your brand of TV is not listed, please use code for the AutoSense. Contact our Technical Support Department at 1-800-328-8908 if none of the codes worked for you.

TV Brand	Three-Digit Device Codes	Four-Digit Device Codes
AutoSense™ (Use this code to allow the Elite Digital Pillow Speaker to automatically identify the TV)	000	0000
RCA® Healthcare TV	001	0001
Zenith® Healthcare TV	002	0002
LG® Healthcare TV	002	0002
PDI® Healthcare TV	003	0003
Philips Healthcare TV	003	0003