



CELEBRATING



YEARS

# Create your successful fall prevention program with Crest Fall Fighter® Monitoring Systems

**CREST**



## Crest Fall Fighter® Family of Fall Management Products

The Fall Fighter family was designed to provide the most critical features for fall prevention in a simple, durable and tamper-proof format. There are three styles of Fall Fighter fall management products: String Monitor, Chair Monitor and Bed Monitor. 1-year and 60-day sensor pads are available for use with the Chair and Bed Monitors. Fall Fighter offers quality features at an affordable price. Choose the monitoring option that fits the needs of your facility.

- “Always On” monitoring
- Multiple alarm tones
- Multiple volume options
- Durable sensor pads with reinforced cords




**CREST**



**TOOLS FOR YOUR SUCCESS!**

Get your free copy of our white paper and Fall Fighter® DVD training video by calling 1-800-328-8908 or visit our website at [www.cresthealthcare.com](http://www.cresthealthcare.com)!

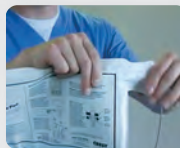

### View Fall Fighter® Videos

Watch the Fall Prevention Case Study Video and Fall Fighter Training Video  
<http://tinyurl.com/26c98pl>



### FREE White Paper!

A case study on one facility's experience making the change to Fall Fighter® products  
<http://tinyurl.com/29keoqk>

Check out what our  
customers have to say  
about Fall Fighter®



“Fall Fighter Monitors make things easier for the staff because there isn't an on/off switch that they have to worry about. The reset is easy to use.”

**Diana, Central Supply**



114448



114687

Part No.	Description
114448	Fall Fighter Chair Monitor
114449	Fall Fighter Bed Monitor
114452	AC Adapter for Fall Fighter Bed Monitor
114687	<b>NEW</b> Monitor Protector for 114448 and 114449 Fall Fighter Monitors
101968	Nurse Call Adapter Cable, 1/4" to 1/4" plug, 7 ft., grey

**GUARANTEED IN-STOCK**  
OR WE'LL SHIP FREIGHT-FREE GROUND

**1**  
YEAR WARRANTY  
ON MONITORS



114685

114479

**Buy as a kit and save more!**

Part No.	Description
114477	Fall Fighter Chair Monitor and 1-Year Chair Sensor Pad
114479	Fall Fighter Bed Monitor and 1-Year Bed Sensor Pad
114685	Fall Fighter Chair Monitor and 60-Day Chair Sensor Pad
114686	Fall Fighter Bed Monitor and 60-Day Bed Sensor Pad

**Crest Fall Fighter® Sensor Pad Monitors**

Fall Fighter® has the robust features you need for chair or bed applications at an affordable price. The alarm automatically resets when patient returns to the pad or is quickly silenced with the reset button. A hinged battery cover with optional screw controls and prevents tampering. The front LED flashes blue to indicate low battery and red to indicate active status and alarm. Durable steel mounting clip included. The Bed Monitor includes a nurse call interface jack, optional AC adapter input and remote-only volume option to silence alarm in room. Uses three AAA batteries (sold separately).

- "Always On" monitoring
- Volume range from 85-105dB
- Chime, warble and buzz alarm tones
- 0-2 second delay reduces false alarms
- Optional monitor protector available to help prevent damage from accidental drops
- AAA batteries have 4 times the current capacity of 9V batteries for longer life



ROBUST PLUG DESIGN

DURABLE CORD CONNECTION FOR LONGER LIFE

NO HOLES FOR BETTER PROTECTION AGAINST INCONTINENCE

**CAUTION!**  
DO NOT BEND, FOLD, OR ROLL SENSOR PAD!

114451

114503

**GUARANTEED IN-STOCK**  
OR WE'LL SHIP FREIGHT-FREE GROUND

**LATEX-FREE**

**1**  
DAY WARRANTY

**GUARANTEED IN-STOCK**  
OR WE'LL SHIP FREIGHT-FREE GROUND

**LATEX-FREE**

**60**  
DAY WARRANTY

**Crest Fall Fighter® Universal Sensor Pads**

Fall Fighter Sensor Pads are built to last and have a 60-day or 1-year warranty. Chair pads are 8" X 14" with corner cord exit to easily direct it through chair back or side. Bed pads are 10" X 30" with side cord exit.

- Over-molded plug end to protect against strain and clip breakage
- Heavy-duty bonding with reinforced anchor prevents separation of cord from pad

Part No.	Description
114503	Fall Fighter Chair Sensor Pad, 60-Day
114504	Fall Fighter Bed Sensor Pad, 60-Day
114450	Fall Fighter Chair Sensor Pad, 1-Year
114451	Fall Fighter Bed Sensor Pad, 1-Year

**Buy Sensor Pads in money-saving 10 packs**

Part No.	Description
114643	Fall Fighter Chair Sensor Pad, 60-Day, 1 pack (10/pk)
114644	Fall Fighter Bed Sensor Pad, 60-Day, 1 pack (10/pk)
114645	Fall Fighter Chair Sensor Pad, 1-Year, 1 pack (10/pk)
114646	Fall Fighter Bed Sensor Pad, 1-Year, 1 pack (10/pk)



LOW BATTERY INDICATOR

BATTERY CHECK BUTTON

HIGH, LOW OR CHIME TONES

114053

**Crest Fall Fighter® String Monitor**

The Fall Fighter® string monitor offers big features including no On/Off switch, multiple volume options and snap-in battery compartment door with optional screw. Check battery power with the test button and low battery LED. Magnet connection strength ensures individuals cannot walk away with the monitor still attached. Includes holster, mounting screws, strap and two AAA batteries.

- "Always On" monitoring
- Battery saving technology
- High, Low or Chime tone
- AAA batteries have 4 times the current capacity of 9V batteries for longer life

Part No.	Description
114053	Fall Fighter String Monitor
113823	Replacement Cord and Magnet
113827	Replacement Holder with strap
114234	Replacement Battery Cover