



Connecting Elite Pillow Speakers to Hospital Grade TVs

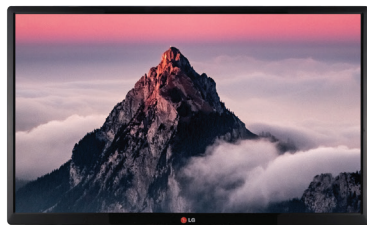

Hospital Grade TV

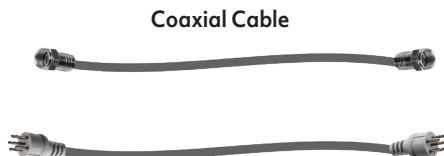

Coaxial Cable

TV Jumper Cable

Plugs into the pillow speaker port of the hospital grade TV



Antenna Coaxial Cable
Cable in-wall,
coaxial jack behind plate



**Crest Elite Digital™
Pillow Speaker
or
Crest Elite Infrared™
Pillow Speaker**



Patient Station

These Crest parts are used to connect a hospital grade TV to a Crest Elite Digital or Infrared pillow speaker. This configuration provides TV audio and remote control function via the pillow speaker along with nurse call and optional light or other environmental controls.