

Birdie™ EVO XPLUS

Invacare® Premier Series Floor Lift



Optimized design for safe and comfortable transfers

SLOW'R™ INTEGRATED DAMPNER | ERGONOMIC HOOK DESIGN



Yes, you can.®

An Evolution in Floor Lifts

The difference is in the details. The new Premier Series **Birdie^{EVO XPLUS}** was designed to offer optimum functionality and an innovative feature set that provides maximum comfort and security when lifting or transferring from a bed, chair or off the floor.

The sleek and modern design was driven by customer insights, with a goal of incorporating what was most valued in a floor lift while bringing in new features to enhance the user and caregiver experience.

Impressive lifting range

The **Birdie^{EVO XPLUS}** offers an impressive lift range that allows for lifting from floors and high surfaces. In addition, the 360° rotation helps make positioning easier for the caregiver.

Easy to assemble, transport and store

The space-saving design allows for toolless assembly, and can be easily folded and unfolded for transport and storage.



Easy to assemble, transport and store

The new **Birdie^{EVO XPLUS}** can be set up, folded and unfolded without the need for tools. It takes up minimal space when stored and is easy to maneuver and transport.



NEW **EXCLUSIVE**
FEATURES

Slow'R®

▶ The boom was designed with an integrated dampener to help reduce rocking, allowing for patient comfort and steady caregiver control during lifts and transfers.



Hook Design

▶ A new ergonomic hook design allows the sling loop to be attached with just one hand. A wider slot makes positioning loops easier, and the hook's nose helps prevent the sling loop from sliding out unintentionally.

Additional NEW features



▶ **Ergonomic push bar** makes handling and maneuvering easy.



▶ **Distinctly designed curvature of the boom and spreader bar** helps caregivers maintain eye contact with patients and residents throughout the transfer.



▶ **Ergonomic leg design** features natural edges and smooth curves, allowing for easy movement around chairs and wheelchairs.




▶ **Premium Casters** maximize maneuverability - especially on challenging surfaces.

Accessory Options ▼

- ▶ Extensive selection of Invacare Sling styles, fabrics and sizes
- ▶ Scale attachment, Item #BLS5
- ▶ Spreader bar padding, Item #AP1660511
- ▶ Protective leg covers, Item #AP1648551
- ▶ Service Data toolkit, Item #AP1541967



Specifications ▼

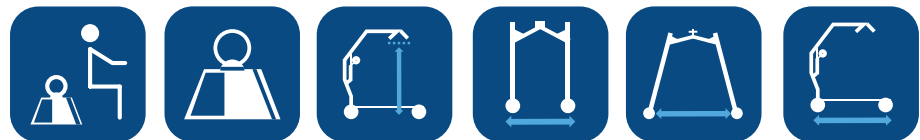

model # ▼	Weight Capacity ▼	Product weight without sling or accessories ▼	Height at Sling Hookup ▼	Base Width (Closed) ▼	Inside Width (Open) ▼	Base length ▼
-----------	-------------------	-----------------------------------------------	--------------------------	-----------------------	-----------------------	---------------

Birdie EVO XPLUS	BRDXPLUS	463 lb.	96.4 lb.	24" - 74" H	25.3"	40.1"	50.4"
-------------------------	----------	---------	----------	-------------	-------	-------	-------



Clearance ▼	Folded Height ▼	Emergency Controls ▼	Battery ▼	Lifts per Charge* ▼	Limited Warranty ▼
-------------	-----------------	----------------------	-----------	---------------------	--------------------

Birdie EVO XPLUS	4.5"	18.5"	Stop, Up and Down	24V DC 2 w/Unit	up to 40 cycles	5 Years Lift 2 Years Electronics
-------------------------	------	-------	-------------------	--------------------	--------------------	-------------------------------------

* With battery charge between 100% to 50%

Regulatory data ▼



Telephone: 800-333-6900
orders@invacare.com
www.pro.invacare.com

Invacare Corporation
 One Invacare Way
 Elyria, OH 44035
 United States



Yes, you can.®