

# PillowCase by HATCH+MED

# CREST HEALTHCARE®

Transform your iPad into an FDA medical device - PillowCase keeps the iPad always powered, secured, and safe from fluid or chemical damage around the patient.




Compatible with iPad 7-10th Gen,  
iPad Air 3-5th Gen  
& iPad mini



Optional kickstand with a  
spring loaded wall & bed rail  
clip



UL listed plug & play  
integration with any nurse call  
system

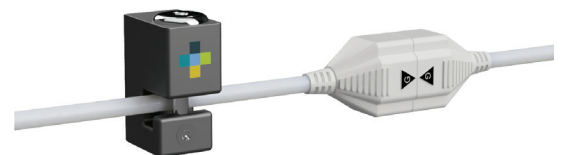

Wall mounted lock allows  
magnet to disconnect while  
preventing iPad theft



info@hatchmed.com

HATCH+MED



844-428-2463

# PillowCase Configurations

## Version A: iPad Security Case



- No physical nurse call buttons on case
- No software integration required
- Direct connection to power outlet with UL listed power supply
- Magnetic cable with 8ft on case side & 2ft on wall side
- Water-resistant and cleaning agent proof enclosure
- Kickstand & security lock accessories

## Version B: Emergency Call Button



- UL 1069 hardware button for emergency calls on any nurse call system via 1/4" port
- No software integration required
- Wall mounted hub with UL listed nurse call connection and power supply
- Magnetic cable with 8ft on case side & 2ft on wall side
- Water-resistant and cleaning agent proof enclosure
- Kickstand & security lock accessories

## **COMING SOON:** Digital Pillow Speaker



- UL 1069 hardware button for emergency calls
- Software room controls, TV controls, clinical requests directly through nurse call system patient station (no API setup required) for Responder 5, HillRom and other systems
- Wall mounted hub with UL listed nurse call connection and power supply
- Magnetic cable with 8ft on case side & 2ft on wall side
- Water-resistant and cleaning agent proof enclosure
- Kickstand & security lock accessories



# PillowCase Specifications

## Size & Weight

	Larger Case	Smaller Case
<b>Supported iPads</b>	iPad 7th, 8th, 9th & 10th Gen iPad Air 3rd, 4th & 5th Gen iPad Pro 10.5" & 11"	iPad mini 6th Gen
<b>Weight</b>	1.0 lb (without iPad) 2.1 lb (with iPad)	0.8 lb (without iPad) 1.5 lb (with iPad)
<b>Dimensions</b>	11.5" W x 7.9" H x 1.3" D	9.7" W x 6.3" H x 1.3" D

## General Info

<b>Connection &amp; Charging</b>	10ft keyed magnetic cable (8ft on PillowCase side & 2ft on wall side)
<b>Materials</b>	Polycarbonate & TPE with antimicrobial & UV stabilizing additives
<b>Fluid Intrusion</b>	IP 65 - dust & water stream resistant
<b>Cleaning</b>	Stable with standard hospital surface cleaners
<b>Power</b>	12VDC @ 3A plug-in power supply provided
<b>iPad Features</b>	Wifi, Bluetooth, cameras, speakers, microphones and screen not affected

## Regulatory

<b>Certifications (Version A)</b>	UL/cUL 62368-1 (power supply) & FDA Registered Medical Device, Class 2 Exempt
<b>Certifications (Versions B &amp; C)</b>	UL/cUL 1069 (nurse call) & FDA Registered Medical Device, Class 2 Exempt
<b>Patents</b>	US10863012B2 US20220066525A1

## Accessories

<b>Kickstand &amp; Rail Clip</b>	Friction hinged kickstand & springed rail clip compatible with all bed & stretcher rails as well as a wall mount
<b>Security Lock</b>	Compatible with all versions, universal key provided

