



## Help Everyone Breathe Easier...

### Nurse Call Ventilator Remote Alarm System - OEM

This fully supervised Nurse Call Ventilator Remote Alarm System by Engineered Electronics provides a simple, reliable method of communicating a ventilator alarm condition through an auxiliary nurse call system input. This supervised ventilator alarm will alert the appropriate staff at the earliest moment if the monitor cable connection is broken at the nurse call station, the vent itself, or anywhere in between, as well as when alarm conditions are generated by the ventilator. Once the ventilator alarm discontinues, the call to the nurse call system will cease. The alarm monitor includes two LED indicators: red for alarm and green for normal operating conditions. The integrated 2-conductor jack allows the continued use of another device that needs connection to the auxiliary input for the nurse call system. The monitor is powered by a TARDIRAN high energy lithium battery (included).

- Plug in connections are supervised outside the nurse call system
- Automatic call when connection is broken
- Unit comes standard with RJ11 connector, additional adapters available
- No mounting required
- 6 second delay cycle




**Safe & Reliable**




Part No.	Description	Price
*900-621	Ventilator Remote Alarm System, supervised	170.95
*500-553	Adapter, DB-9 plug, for Puritan Bennet 7200	42.95
*500-554	Adapter, 1/4" Stereo Plug, for Resprionics and Phillips V60	37.95
*500-555	Adapter, Hirose Plug, for Puritan Bennet 760 and 840	116.95
*500-562	Dummy Plug, RJ11	6.95
*200-007	Replacement Lithium Battery, AA	5.95

\*Item is Drop Shipped. Additional Adapters Available. Please contact Crest for details.



-  800.328.8908
-  800.369.9207
-  CrestHealthcare.com

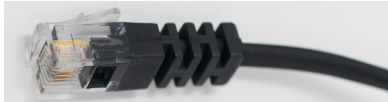

©2022 Crest Healthcare. \*Shipped direct from the supplier to you. All listed prices and products are subject to change without notice. Crest is not responsible for photographic or typographical errors.