



■ Nexus Master Station

NEXUS

Nexus is a **Wireless Call System** based around an Android display unit. It has been designed from ground up to be **Flexible, Easy To Program, Intuitive To Use & Competitively Priced.** It is great for any small facility with up to 25 beds. Although the system is designed for ease-of-use, Nexus is a professionally designed system with full functionality. Nexus is developed based on UL 2560. The system is fully supervised to provide a reliable solution.



Touchscreen Master Station

- Color Coded Touchscreen Display
- Calls Ordered By Priority
- Audio Call Levels Indication
- Low Battery And Device Lost Alerts
- Event Log and Response Time Tracking
- Password Protected
- Options for Wall-mount / Desk-mount

NEXUS Key Features & Benefits :

-  **DIY Nurse Call System**
Plug & Play structure, easy to program, simple installation
-  **Easy-To-Use**
Intuitive UI design, straightforward to user and no specialized training required
-  **Easy Deployment**
Fully functional nurse call system can be deployed within a matter of minutes
-  **OverDoor Indicator**
Visual and audible alarm, when acts as a signal repeater for extending coverage
-  **Fault Supervision**
Alerts for product faults, out of range and low battery conditions
-  **Event Log**
Every event that occurs on the system is collected. Data can be exported for analysis and reporting
-  **Reliable Backups & Rapid Recovery Capabilities**
Instant recovery solution in the event of an unexpected hardware or system failure
-  **Messaging Solutions Available**
Instant notification sent to staff members via in-house pagers or smartphones
-  **Real-time Remote Display Panels**
Staff do not have to travel far to see the location of the call
-  **Fall / Wander Management Integration**
Integrates seamlessly with the fall / wander management solutions. Direct alerts to the nurse call system for central monitoring and response. Local alerts can be disabled for silent alarming

Assisted Living Facility



Fall Management



Fall Monitor

The CordFree monitor acts as a room controller unit. Manages up to 9 sensors and supports a wide range of sensor devices



Sensor Pad

Various types & sizes of sensor pads, include above or under-mattress bed pads, chair pads, floor mats are available



Portable Handset

A portable nurse call handset can be fastened on a bed rail. It provides the most affordable solution for bedding needs




Bedhead Call Point

Call & Reset 2-in-1 bedhead unit is completely wireless. Battery operated with long battery life expectancy. 4P4C jack for a nurse call handset



Nurse Station

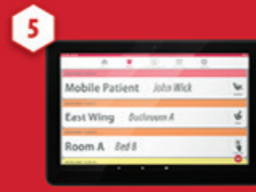

Master Station

A touchscreen Master Station indicates calls in order of priority. Can either stand on a desk or wall-mounted



Wearable Device




Pendant Transmitter

The Pendant Transmitter can be worn as a wristband or pendant, and allow patients to call for emergency alarms. Automatic fall-alert is available for pendant use



Wander Management

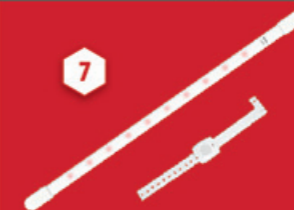

Door Strip Sensor & Wristband Transmitter

Provides audible and visual alert when patient wearing wristband transmitter move within the range of protected doors



Bathroom / WC Safety



OverDoor Indicator

Provides visual and audible alarm, while acting as a signal repeater for extending coverage



Bathroom Call Point






Wireless Call & Reset 2-in-1 button is waterproof to IP65. It comes with a pull cord to ensure accessibility from the toilet or floor

Multiple Options • Proven Performance

System Components

 <p>NDP-21</p>	<p>Master Station 10.1" Rugged Tablet PC • With Nexus Software Installed • IP65 • 8000mAh Battery</p>	 <p>ILB-21</p>	<p>Overdoor Indicator Dual Colored Indicator Light • Buzzer With Adjustable Volume • Signal Repeater • Battery Back-up In Case Of Power Low</p>
--	--	--	--

Nurse Call Components

 <p>WCP-11C</p>	<p>Bedhead Call Point Call & Reset Buttons 2-in-1 • Call Reassurance LED • 4P4C Connection For Nurse Call Handset (NEC-12) • Replaceable Batteries • Long Battery Life</p>	 <p>NEC-12</p>	<p>Nurse Call Handset Easy To Press Big Red Call Button • Call Reassurance LED • Standard Cable With USB Break Point & 4P4C Connector • Comes With Wall-mounted Bracket • Compatible With WCP-11C</p>
 <p>WCP-11</p>	<p>Bathroom / WC Call Point Call & Reset Buttons 2-in-1 • Comes With Removable Pull Cord • IP65 Waterproof Rating • Replaceable Batteries • Long Battery Life</p>	 <p>NEC-05</p>	<p>Portable Handset Easy To Press Big Red Call Button • Green Button To Reset • Call Reassurance LED • 2xAAA Batteries Operated • Comes With Clip & Plastic Wall-mounted Bracket for Storage</p>
 <p>TXP-11</p>	<p>Pendent Transmitter Active Emergency Button • Man-down Detection For Pendant • Waterproof to IP65 • Replaceable Batteries • Includes Lanyard With Anti-ligature Point • Optional Watch Strap (Not Included)</p>		

Fall Management Components

 <p>NGM-21</p>	<p>CordFree Monitor Monitor Up to 9 Sensor Devices • Pause Mode • Sleep Mode • Low Battery Alert • Signal Lost Alert • Tamper-Proof Reset • Delay Setting • 3 Tone Selections • Adjustable Volume • Nurse Call Output • Operates On AC Adaptor Or Batteries • Protective Cover • Hanger Strap</p>	 <p>WPAS-10</p>	<p>Early-Alert Sensor Bed Pad Sends Alarm BEFORE Patient Has Left The Bed • Multi-Zone Detection Design • 30"x10" (76 x 25 cm) • Water Resistant • Built-in Transmitter</p>
 <p>WMBP-20</p>	<p>Cordless Under-Mattress Sensor Pad Works Beneath Most Mattresses • 30" x 10" (76 x 25cm) • Water Resistant • Anti-Microbial Surface • Works With Rondish Transmitter (TM-11v3)</p>	 <p>CBP-01</p>	<p>Cordless Bed Pad 30" x 20" (76 x 51cm) • Water Resistant • Works With Rondish Transmitter (TM-11v3)</p>
 <p>CMAT-02</p>	<p>Cordless Floor Mat 48" x 24" (122 x 61cm) • Heavy Duty • Non-Slip Material • Works With Rondish Transmitter (TM-11v3)</p>	 <p>CCP-01</p>	<p>Cordless Chair Pad 15" x 10" (38 x 25cm) • Water Resistant • Works With Rondish Transmitter (TM-11v3)</p>
 <p>CFMT-10</p>	<p>Cordless Foldable Floor Mat 36" x 24" (90 x 61cm) • Foldable Design • Heavy Duty • Non-Slip Material • Works With Rondish Transmitter (TM-11v3)</p>	 <p>PIR-21</p>	<p>Passive Infrared Sensor LED Status Indicator • Internal Alarm Or Signal Transmit Options • Pause Mode • Easy Reset Button • Adjustable Detection Range • 3xAA Battery Operated • Low Battery Alert • DC (6V) Input Option</p>

Wander Management Components

 <p>DMS-21</p>	<p>Door Strip Sensor Visual & Audible Alert • Adjustable Sensing Range • Automatic Door Lock Option • Customize ID • Fault-checked</p>	 <p>TXP-21</p>	<p>Wristband Transmitter Unique code to identify each tag • Tag can be switched off using magnetic switch • Waterproof to IP65 • 3-year battery life with low battery alert • Comes with faux leather watch strap</p>
--	---	--	--

Optional Components

 <p>DOCK-01</p>	<p>Docking Station 3 USB Extend Docking • Support Ethernet</p>
---	---



We aim to improve the quality of life for elderly, patients and caregivers with human-centric technology.



Tel. : (852) 2541 6685
 Fax : (852) 2541 7411
 Email : sales@rondish.com

Website : www.rondish.com
 Facebook : Facebook/RondishCare